



## Easi-Scan Curve

Applications



### Curved for a deeper and wider image

Easi-Scan Curve is ideal for ageing pregnancies around 60 days and until the fetus drops over the pelvic brim (about 120 days). It also has numerous non-reproductive applications including abdominal examinations of other ruminant species and small animals.

#### Better calving date estimation in late gestation

The Easi-Scan Curve probe penetrates much deeper than the traditional linear probe. This allows you to visualise later pregnancies and determine fetal age by measuring eye socket diameter up to 8 months. This promotes easier farm management and better herd health.

#### Excellent image quality

The probe has a 60 mm radius and 60 percent curve to give you a wide field of view and up to 24 cm depth, with detailed image information at near and deep range. The 128 elements probe produces 40 frames per second across 16 digital channels, so there is no lag or image loss when scanning rapidly.

#### Comfortable working

Easi-Scan is light. You carry it comfortably on your hip. Your hands are free and there is no strain on your neck, shoulders or back. When used with the BUG goggle it allows great peripheral vision.

Visit our website for product in use video, full technical specifications and clinical learning materials.

Four buttons. Easy to operate.

Probe is fixed and sealed for durability and less downtime caused by connection failures and dust/water penetration.

Probe designed from the ground up, curved for greater depth (8, 12, 16, 24 cm mode)


Detachable battery. Lasts up to 7 hours



37 day pregnancy  
8 cm mode \*



94 day calf head  
24 cm mode \*

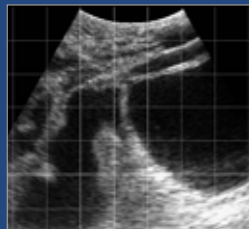

Cotyledons  
24 cm mode \*



Uterine horns, CL  
12 cm mode \*

easi-scan™  
curve

\* Zoomed for printing

### BCF because ...

We help you grow your skills with practical on-farm courses and the largest online learning centre.

You can rely on BCF. We've been in animal imaging for over 30 years, our products are made in Scotland, UK and supported by a strong network of more than 40 specialist veterinary distributors worldwide.



### Viewing options and accessories



BUG Goggles



Smart Display



Remote Display



Carry Case



Wrist monitor



Easi-Scan Grip



Introducer



Gestational age app