

Midmark Professional Sterilization Systems



60 Litre B Class Sterilizer

Version : 02/2013



Standard Features


- High throughput—large 60 litre capacity
- Built in printer and USB connection
- Class B with EN 13060 compliance
- Wide range of cycles at 121°C and 134°C for porous loads and pouches
- Integrated reverse osmosis system
- Supplied with tray holder and 4 trays
- Low water consumption—3.5 litres per cycle
- Power connection—16A 240V 50Hz-60Hz
- Single phase plug
- User friendly operation

Specifications

- Height: 820mm (1600mm with stand)
- Width: 646mm (680mm with stand)
- Depth: 783mm
- Weight: 160kg (205kg with stand)
- Chamber volume: 59.9 litres
- Chamber diameter: 355mm
- Chamber depth: 605mm
- Free sterilisable weight: 15kg
- Sterilisable package weight: 7kg


STERILIZATION CYCLES

134°C classe B

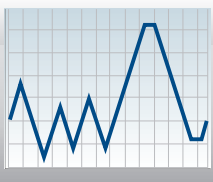


■ Sterilization temperature.....	134°C
■ Sterilization pressure.....	2,1 bar
■ Sterilization time.....	5'
■ Drying time.....	15'
■ Cycle time*.....	1h15'

121°C classe B

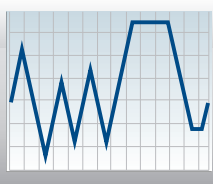


■ Sterilization temperature.....	121°C
■ Sterilization pressure.....	1,1 bar
■ Sterilization time.....	18'
■ Drying time.....	15'
■ Cycle time*.....	1h25'

134°C PRION

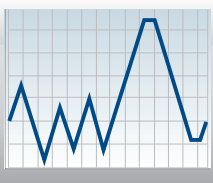

■ Sterilization temperature.....	134°C
■ Sterilization pressure.....	2,1 bar
■ Sterilization time.....	20'
■ Drying time.....	15'
■ Cycle time*.....	1h35'

134°C classe S

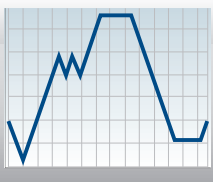

■ Sterilization temperature.....	134°C
■ Sterilization pressure.....	2,1 bar
■ Sterilization time.....	5'
■ Drying time.....	15'
■ Cycle time*.....	1h



* Approx. cycle length calculated with warm chamber at full load.