




## Nasal cannulae

A range of cannulae for adults, paediatrics and neonates

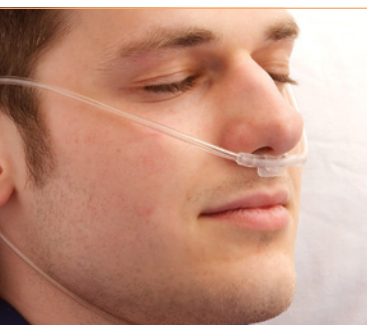

Oxygen and Aerosol Therapy • Nasal Cannulae

## Nasal cannula

### Conventional nasal cannula

A range of products for adult, paediatric and neonates are available with a wide selection of prong shapes to suit all patient requirements. All oxygen tubing in the range is manufactured in a 'multi-channel' design, which permits oxygen flow even if the tubing is kinked.

### Also available in Satin


Cannula made from Intersurgical's soft satin material, proven in testing to create 20% less friction<sup>1</sup> on the patient's face and ears compared to a standard cannulae.

### Prong shape


The following prong shapes are available:




▶ Video available at [www.intersurgical.com](http://www.intersurgical.com)

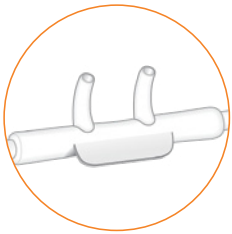


**Curved prongs:**  
Improved anatomical shape and fit

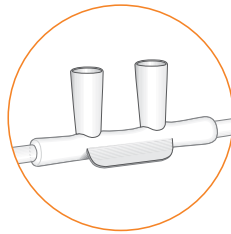

**Flared prongs:**  
Slow down the flow of oxygen

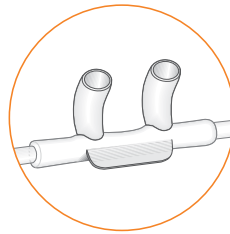

**Curved and Flared prongs:**  
Combines the benefits of both

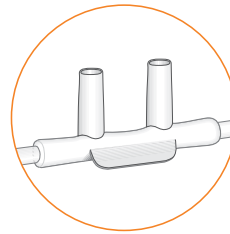

**Straight Prongs**



**Curved prongs with ear guards:**  
Improved anatomical fit and patient comfort

### Now available with ear guards

The Earguard option includes soft silicone foam cushions, which are easily adjusted and securely fitted to prevent trauma to the patient's ears.



### Phthalate free nasal cannula

The Phthalate free range of nasal cannula has been designed to address the concerns of neonatal exposure to DEHP in premature infants and newborns requiring intensive care in a neonatal unit.





**Multi-Channel oxygen tubing**

Allows oxygen to continue to flow if the tubing becomes kinked

**Non-slip toggle**

Enables user to adjust for best possible fit

**Latex free**

Suitable for patients with latex allergy

**Clear material**

Unobtrusive

(Low friction Satin material also available)

**Over-ear design**

Secure positioning

(Ear guards available for added comfort)

**Soft prongs**


Minimise irritation

**Location saddle**


Assists in the optimum positioning of the prongs

**Adult headset version**

Ideal for domiciliary use

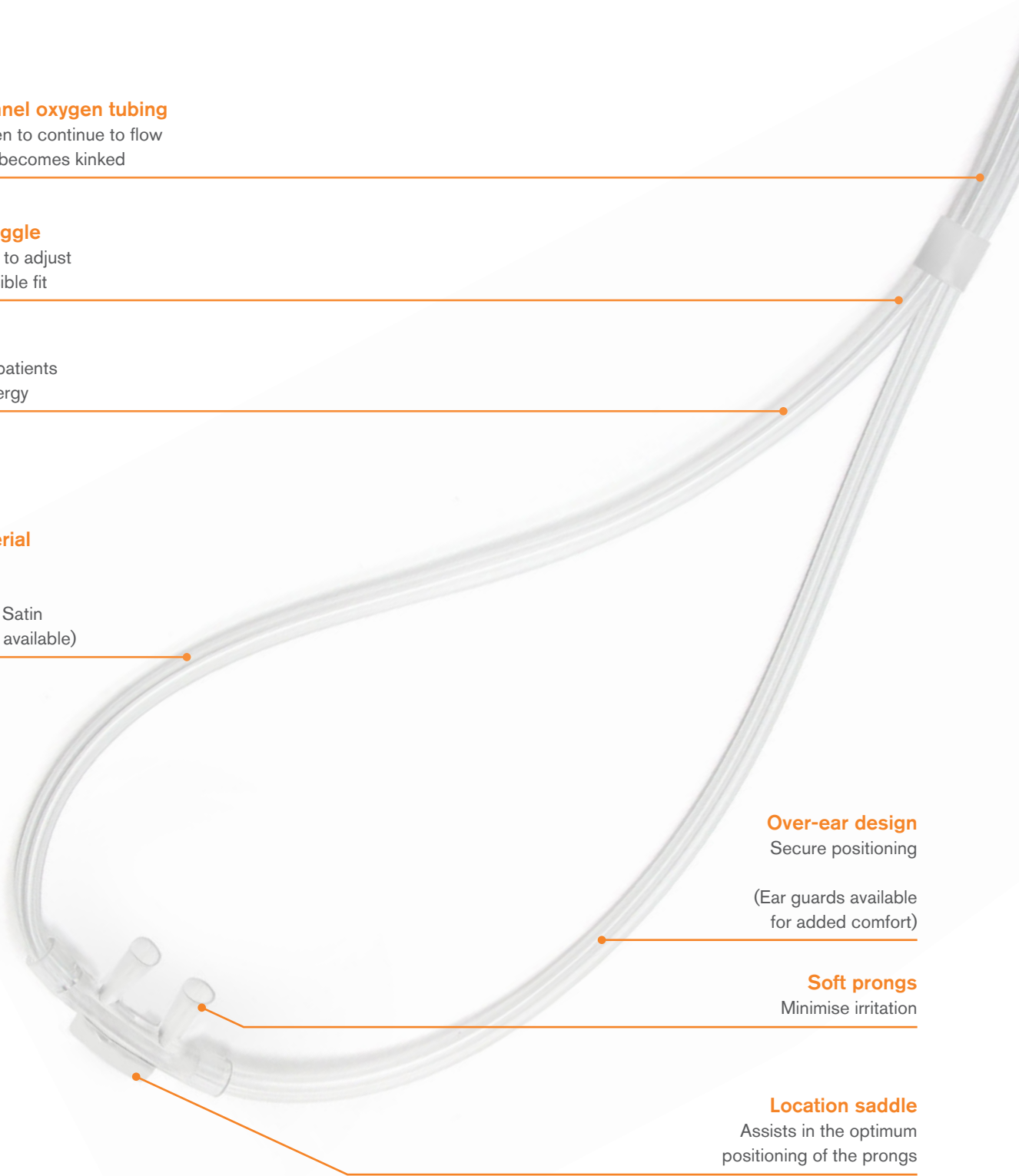

### Adult range of nasal cannula

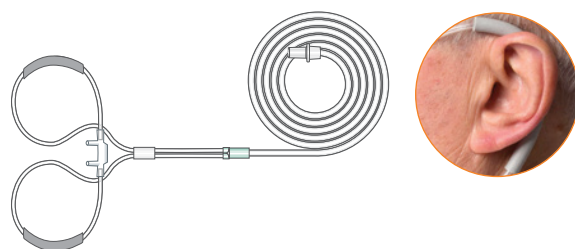

1165002

Code	Description	Length	Box Qty.
1161000	Adult, nasal cannula with straight prongs and tube	1.8m	50
1161004	Adult, nasal cannula with straight prongs and tube	2.1m	50
1162000	Adult, nasal cannula with straight prongs and tube	5m	20
1169000	Adult, nasal cannula with straight prong, headset and tube	0.5m	100
1165000	Adult, nasal cannula with curved prongs and tube	1.8m	50
1161020	Adult, Satin nasal cannula with straight prongs and tube	1.8m	50
1165020	Adult, Satin nasal cannula with curved prongs and tube	1.8m	50
1165002	Adult, nasal cannula with curved prongs, ear guards and tube	2.1m	50
1165003	Adult, nasal cannula with curved prongs, ear guards and tube	5m	40
1165021	Adult, Satin nasal cannula with curved prongs and tube	2.1m	50
1161021	Adult, Satin nasal cannula with straight prongs and tube	2.1m	50
1166000	Adult, nasal cannula with flared prongs and tube	1.8m	50
1167000	Adult, nasal cannula with curved/flared prongs and tube	1.8m	50
1168000	Adult, nasal cannula with curved/flared prongs headset		100

### Neonatal and paediatric range of nasal cannula

Code	Description	Length	Box Qty.
1163000	Paediatric, nasal cannula with curved prongs and tube	2.1m	50
1164000	Neonatal, nasal cannula with curved prongs and tube	2.1m	50

### Phthalate free range of nasal cannula

Code	Description	Length	Box Qty.
1160002	Infant, nasal cannula with curved prongs and tube	2.1m	50
1160001	Neonatal, nasal cannula with curved prongs and tube	2.1m	50
1160000	Premature, nasal cannula with curved prongs and tube	2.1m	50

Phthalate free paediatric, infant and neonatal options



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK  
 T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015

**Please think before you print**  
 Save energy and paper.  
 If you must print this information sheet please print it double sided.